

# Go Direct<sup>®</sup> Gas Pressure

(Order Code GDX-GP)



Go Direct Gas Pressure is used to monitor pressure changes in gas-law experiments in chemistry, biology, and physics.

**Note:** Vernier products are designed for educational use. Our products are not designed nor are they recommended for any industrial, medical, or commercial process such as life support, patient diagnosis, control of a manufacturing process, or industrial testing of any kind.

## What's Included

- Go Direct Gas Pressure
- Two tapered valve connectors inserted into a No. 5 stopper
- One tapered valve connector inserted into a No. 1 stopper
- One two-way valve
- Two Luer-lock connectors connected to either end of a piece of plastic tubing
- One 20 mL syringe
- Two transpiration tubing clamps
- Micro USB Cable

## Compatible Software

See [www.vernier.com/manuals/gdx-gp](http://www.vernier.com/manuals/gdx-gp) for a list of software compatible with Go Direct Gas Pressure.

## Getting Started

Please see the following link for platform-specific connection information:

[www.vernier.com/start/gdx-gp](http://www.vernier.com/start/gdx-gp)

### Bluetooth Connection

1. Install Vernier Graphical Analysis™ on your computer, Chromebook™, or mobile device. If using LabQuest®, make sure LabQuest App is up to date. See [www.vernier.com/ga4](http://www.vernier.com/ga4) for Graphical Analysis availability or [www.vernier.com/downloads](http://www.vernier.com/downloads) to update LabQuest App.
2. Charge your sensor for at least 2 hours before first use.
3. Turn on your sensor by pressing the power button once. The Bluetooth®

### USB Connection

1. Install Graphical Analysis on your computer or Chromebook. If using LabQuest, make sure LabQuest App is up to date. See [www.vernier.com/ga4](http://www.vernier.com/ga4) for software availability or [www.vernier.com/downloads](http://www.vernier.com/downloads) to update LabQuest App.
2. Connect the sensor to the USB port.
3. Launch Graphical Analysis or turn on LabQuest. You are now ready to collect data.

4. LED will blink red.
5. Launch Graphical Analysis or turn on LabQuest.
6. If using Graphical Analysis, click or tap Sensor Data Collection. If using LabQuest, choose Wireless Device Setup > Go Direct from the Sensors menu.
7. Select your Go Direct sensor from the list of Discovered Wireless Devices. Your sensor's ID is located near the barcode on the sensor. The Bluetooth LED will blink green when it is successfully connected.
8. Click or tap Done to enter data-collection mode.

**Note:** This sensor does not work with the original LabQuest. It works with LabQuest 2 or LabQuest 3.

## Charging the Sensor

Connect Go Direct Gas Pressure to the included Micro USB Cable and any USB device for two hours.

You can also charge up to eight Go Direct Gas Pressure Sensors using our Go Direct Charge Station, sold separately (order code: GDX-CRG). An LED on each Go Direct Gas Pressure indicates charging status.

Charging	Orange LED next to the battery icon is solid while the sensor is charging.
Fully charged	Green LED next to the battery icon is solid when the sensor is fully charged.

## Powering the Sensor

Turning on the sensor	Press button once. Red LED indicator flashes when unit is on.
Putting the sensor in sleep mode	Press and hold button for more than three seconds to put into sleep mode. Red LED indicator stops flashing when sleeping.

## Connecting the Sensor

See the following link for up-to-date connection information:

[www.vernier.com/start/gdx-gp](http://www.vernier.com/start/gdx-gp)

### Connecting via Bluetooth

Ready to connect	Red LED next to the Bluetooth icon flashes when sensor is awake and ready to connect.
Connected	Green LED next to the Bluetooth icon flashes when sensor is connected via Bluetooth.

### Connecting via USB

Connected and charging	Orange LED next to the battery icon is solid when the sensor is connected to Graphical Analysis via USB and the unit is charging. LED next to Bluetooth icon is off.
Connected, fully charged	Green LED next to the battery icon is solid when the sensor is connected to Graphical Analysis via USB and fully charged. LED next to Bluetooth icon is off.
Charging via USB, connected via Bluetooth	Orange LED next to the battery icon is solid when the sensor is charging. Green LED next to the Bluetooth icon flashes.

### Identifying the Sensor

When two or more sensors are connected, the sensors can be identified by tapping or clicking Identify in Sensor Information.

### Using the Product

Connect the sensor following the steps in the Getting Started section of this user manual.

**Important:** The Go Direct Gas Pressure sensing element will be damaged with direct contact to liquid.

### Videos

View videos related to this product at [www.vernier.com/gdx-gp](http://www.vernier.com/gdx-gp)

### Calibrating the Sensor

You do not have to perform a new calibration when using Go Direct Gas Pressure. The sensor is calibrated prior to shipping.

If you do choose to calibrate, a one-point calibration at atmospheric pressure is adequate for most applications. Go Direct Gas Pressure has been calibrated to read station pressure. Station pressure is the true atmospheric pressure at your location, or station. If you prefer it to read sea-level pressure for conducting weather studies, you can perform a one-point calibration to correct for elevation, using the sea-level corrected pressure value for your location obtained from a reputable source (e.g., NOAA, Weather Underground, etc.). Sea-level pressure is

the pressure after the station pressure has been adjusted to its equivalent pressure at sea level. This is commonly done to normalize pressures at various altitudes for weather forecasts.

To calibrate the sensor in Graphical Analysis, complete the following steps.

1. Click or tap the sensor meter to view sensor options.
2. Select Calibrate and follow the applicable prompts on the Calibrate Sensor screen.

Once you have calibrated a Go Direct sensor, the calibration is automatically stored to the sensor and will be used each time you connect to your device. You can always choose to restore factory defaults if you feel the custom calibration is invalid.

### Specifications

Pressure range	0 kPa to 400 kPa
Maximum pressure without permanent damage	410 kPa
Typical accuracy	±3 kPa
Internal volume	0.8 mL
Resolution	0.03 kPa
Maximum sampling rate	50 samples/s
USB specification	2.0
Wireless specification	Bluetooth 4.2
Maximum wireless range	30 m
Dimensions	9.5 cm × 6.0 cm × 3.25 cm
Battery	300 mA Li-Poly
Battery life (single full charge)	~24 hours with active use
Battery life (long term)	~500 full charge cycles (several years depending on usage)

### Care and Maintenance

#### Battery Information

Go Direct Gas Pressure contains a small lithium-ion battery. The system is designed to consume very little power and not put heavy demands on the battery. Although the battery is warranted for one year, the expected battery life should be several years. Replacement batteries are available from Vernier (order code: GDX-BAT-300).

## Storage and Maintenance

To store Go Direct Gas Pressure for extended periods of time, put the device in sleep mode by holding the Power button down for at least three seconds. The red LED will stop flashing to show that the unit is in sleep mode. Over several months, the battery will discharge but will not be damaged. After such storage, charge the device for a few hours, and the unit will be ready to go.

Exposing the battery to temperatures over 35°C (95°F) will reduce its lifespan. If possible, store the device in an area that is not exposed to temperature extremes.

## Water Resistance

Go Direct Gas Pressure is not water resistant and should never be immersed in water.

If water gets into the device, immediately power the unit down (press and hold the power button for more than three seconds). Disconnect the sensor and charging cable, and remove the battery. Allow the device to dry thoroughly before attempting to use the device again. Do not attempt to dry using an external heat source.

## How the Sensor Works

The sensor in this unit has a membrane that flexes as pressure changes. This sensor is arranged to measure absolute pressure. One side of the membrane is a vacuum, while the other side is open to the atmosphere. The sensor produces an output voltage that varies in a linear way with absolute pressure. It includes special circuitry to minimize errors caused by changes in temperature.

## Troubleshooting

For troubleshooting and FAQs, see [www.vernier.com/till/3850](http://www.vernier.com/till/3850)

## Repair Information

If you have watched the related product video(s), followed the troubleshooting steps, and are still having trouble with your Go Direct Gas Pressure, contact Vernier Technical Support at [support@vernier.com](mailto:support@vernier.com) or call 888-837-6437. Support specialists will work with you to determine if the unit needs to be sent in for repair. At that time, a Return Merchandise Authorization (RMA) number will be issued and instructions will be communicated on how to return the unit for repair.

## Accessories/Replacements

### Item

Gas Pressure Sensor Replacement Parts

#1 1-hole Rubber Stopper

#5 2-hole Rubber Stopper

Luer-Lock Connector

Plastic 2-Way Valve

Plastic Syringe

### Order Code

PS-ACC

PS-STOP1

PS-STOP5

PS-LUER

PS-2Way

PS-SYR

### Item

Plastic Tubing

Plastic Tubing Clamps (Pack of 100)

Stopper Stem

Vernier Micro USB Cable

USB-C to Micro USB Cable

Go Direct 300 mAh Replacement Battery

### Order Code

PS-TUBING

PTC

PS-STEM

CB-USB-MICRO

CB-USB-C-MICRO

GDX-BAT-300

## Warranty

Warranty information for this product can be found on the Support tab at [www.vernier.com/gdx-gp](http://www.vernier.com/gdx-gp)

General warranty information can be found at [www.vernier.com/warranty](http://www.vernier.com/warranty)

## Disposal

When disposing of this electronic product, do not treat it as household waste. Its disposal is subject to regulations that vary by country and region. This item should be given to an applicable collection point for the recycling of electrical and electronic equipment. By ensuring that this product is disposed of correctly, you help prevent potential negative consequences on human health or on the environment. The recycling of materials will help to conserve natural resources. For more detailed information about recycling this product, contact your local city office or your disposal service.

Battery recycling information is available at [www.call2recycle.org](http://www.call2recycle.org)

Do not puncture or expose the battery to excessive heat or flame.



The symbol, shown here, indicates that this product must not be disposed of in a standard waste container.

## Federal Communication Commission Interference Statement

This equipment has been tested and found to comply with the limits for a Class B digital device, pursuant to Part 15 of the FCC rules. These limits are designed to provide reasonable protection against harmful interference in a residential installation. This equipment generates, uses and can radiate radio frequency energy and, if not installed and used in accordance with the instructions, may cause harmful interference to radio communications. However, there is no guarantee that interference will not occur in a particular installation. If this equipment does cause harmful interference to radio or television reception, which can be determined by turning the equipment off and on, the user is encouraged to try to correct the interference by one or more of the following measures:

Reorient or relocate the receiving antenna.

Increase the separation between the equipment and receiver.

Connect the equipment into an outlet on a circuit different from that to which the receiver is connected.

Consult the dealer or an experienced radio/TV technician for help.

### FCC Caution

This device complies with Part 15 of the FCC Rules. Operation is subject to the following two conditions:

(1) this device may not cause harmful interference and

(2) this device must accept any interference received, including interference that may cause undesired operation

### RF Exposure Warning

The equipment complies with RF exposure limits set forth for an uncontrolled environment. The antenna(s) used for this transmitter must not be co-located or operating in conjunction with any other antenna or transmitter. You are cautioned that changes or modifications not expressly approved by the party responsible for compliance could void your authority to operate the equipment.

## IC Statement

This device complies with Industry Canada license-exempt RSS standard(s). Operation is subject to the following two conditions:

(1) this device may not cause interference, and

(2) this device must accept any interference, including interference that may cause undesired operation of the device.

**Industry Canada - Class B** This digital apparatus does not exceed the Class B limits for radio noise emissions from digital apparatus as set out in the interference-causing equipment standard entitled "Digital Apparatus," ICES-003 of Industry Canada. Operation is subject to the following two conditions: (1) this device may not cause interference, and

(2) this device must accept any interference, including interference that may cause undesired operation of the device.

To reduce potential radio interference to other users, the antenna type and its gain should be so chosen that the equivalent isotropically radiated power (e.i.r.p.) is not more than that permitted for successful communication.

**RF exposure warning:** The equipment complies with RF exposure limits set forth for an uncontrolled environment. The antenna(s) used for this transmitter must not be co-located or operating in conjunction with any other antenna or transmitter.

Le présent appareil est conforme aux CNR d'Industrie Canada applicables aux appareils radio exempts de licence. L'exploitation est autorisée aux deux conditions suivantes :

(1) l'appareil ne doit pas produire de brouillage, et

(2) l'appareil doit accepter toute interférence radioélectrique, même si cela résulte à un brouillage susceptible d'en compromettre le fonctionnement.

Cet appareil numérique respecte les limites de bruits radioélectriques applicables aux appareils numériques de Classe B prescrites dans la norme sur le matériel interférant-brouilleur: "Appareils Numériques," NMB-003 édictée par Industrie Canada. L'utilisation est soumise aux deux conditions suivantes:


(1) cet appareil ne peut causer d'interférences, et

(2) cet appareil doit accepter toutes interférences, y comprises celles susceptibles de provoquer un dysfonctionnement du dispositif.

Afin de réduire les interférences radio potentielles pour les autres utilisateurs, le type d'antenne et son gain doivent être choisis de telle façon que l'équivalent de puissance isotrope émise (e.i.r.p.) n'est pas plus grand que celui permis pour une communication établie.

**Avertissement d'exposition RF:** L'équipement est conforme aux limites d'exposition aux RF établies pour un environnement non supervisé. L'antenne (s) utilisée pour ce transmetteur ne doit pas être jumelée ou fonctionner en conjonction avec toute autre antenne ou transmetteur.

**Note:** This product is a sensitive measurement device. For best results, use the cables that were provided. Keep the device away from electromagnetic noise sources, such as microwaves, monitors, electric motors, and appliances.

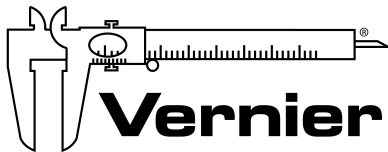

Vernier Software & Technology  
13979 SW Millikan Way • Beaverton, OR 97005-2886  
Toll Free (888) 837-6437 • (503) 277-2299 • Fax (503) 277-2440  
info@vernier.com • www.vernier.com

Rev. 2/18/2021

Go Direct, Vernier Graphical Analysis, LabQuest, and other marks shown are our trademarks or registered trademarks in the United States. All other marks not owned by us that appear herein are the property of their respective owners, who may or may not be affiliated with, connected to, or sponsored by us.

The Bluetooth® word mark and logos are registered trademarks owned by the Bluetooth SIG, Inc. and any use of such marks by Vernier Software & Technology is under license. Other trademarks and trade names are those of their respective owners.